

Anigen Rapid FCoV Ag

General Information

Purpose : Detection of Feline Corona virus antigen

Principle : Immunochromatographic assay

Specimen : Feline feces

Materials provided :

- Anigen Rapid FCoV test Device
- Specimen tubes containing assay diluent buffer
- Sample collection swabs
- Disposable droppers

Special Features :


- One-step testing procedure : Quickly & Accurately
- High sensitivity and specificity
- Sensitivity : 95.5 % vs. PCR
- Specificity : 98.8 % vs. PCR



Feline infectious peritonitis ; FIP

There are two basic forms of FIP-effusive(wet) and FIP-noneffusive(dry) have been characterized. Cats with effusive FIP usually only survive a few days-weeks at best. Cats with noneffusive FIP can survive for weeks or months, although after neurologic signs begin death usually occurs rapidly.

Test Procedure

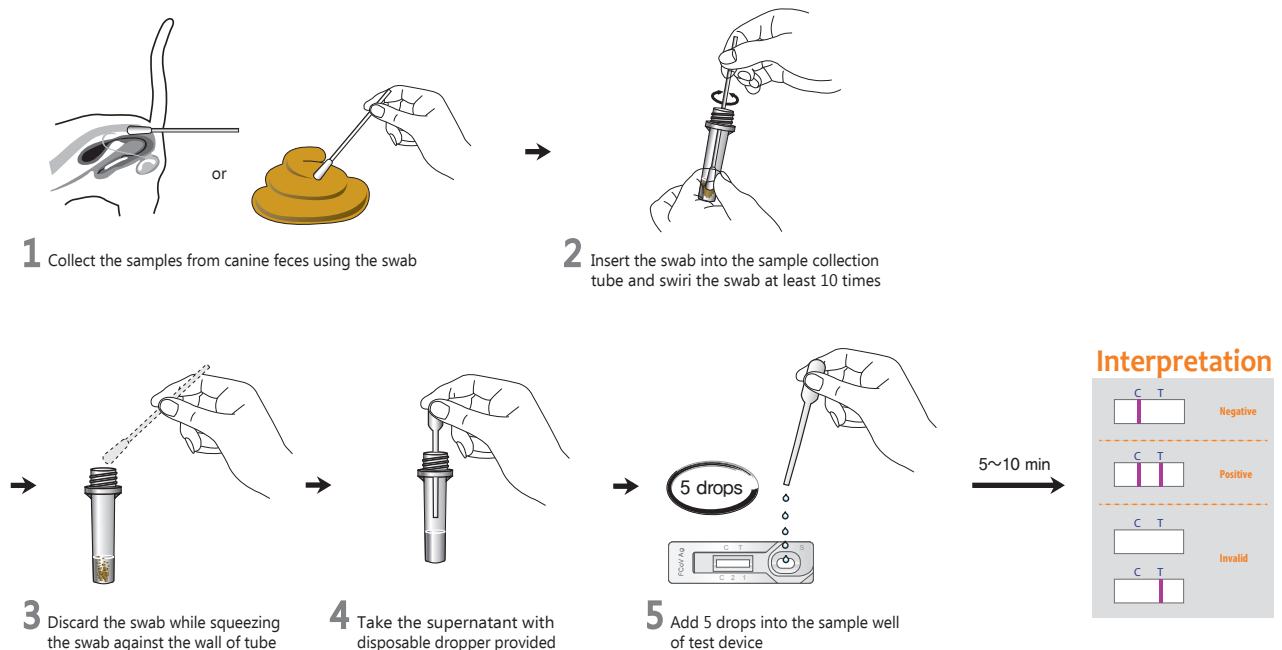

Ordering Information

Cat. No.	Description	Type	Pack size
RG12-05	FCoV Ag	Device	10 Tests/Kit